

DNA
DIAGNOSTIC

Salmonella Velox

Salmonella Velox

Can do it faster...

Why wait for results?

- Results in only 5.5 - 9 hours

Foods in 5.5 hours (incl. enrichment) and primary production samples in 9 hours (Incl. enrichment)

- Nordval International certified

Beef, pork, poultry, fish & seafood, primary production samples and environmental testing

- Species differentiation

Choose to detect all Salmonella spp. and/or S. Enteritidis and S. Typhimurium

- Cost efficient

High throughput and fast product release

- High sensitivity

Detects 1 CFU in 25 g sample

Salmonella Velox

Fastest qPCR test
on the market

Introduction

All Salmonella Velox and Salmonella Velox SE+ST* products are NordVal certified as rapid Salmonella qPCR tests for testing the following categories: Raw and processed meat, poultry and fish/seafood, food swabs, primary production samples and environmental samples.

*Salmonella Enteritidis (SE) + Salmonella Typhimurium (ST)

Advantages

The Salmonella Velox product line has a shorter enrichment time, and less enrichment media, compared to other conventional Salmonella tests.

Salmonella Velox is the fastest certified salmonella test on the market.

Salmonella Velox products run in 96 well formats. All steps are optimized for easy handling. Velox products are easy to use in high-throughput screenings.

Salmonella Velox products provides time and cost-saving systems for Salmonella spp. detection and/or S. Enteritidis and S. Typhimurium.


Total Solutions

All DNA Diagnostic kits are supplied with all necessary reagents and consumables, including:

- Enrichment broth
- Plastic consumable incl. Enrichment bags
- DNA extraction kits
- Pre-aliquoted qPCR mastermix

Description

SV960SEST

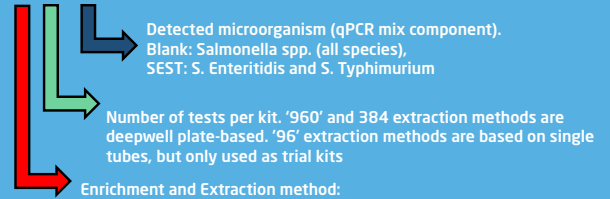

SV: Salmonella Velox (standard), the method for Food samples

SP: Salmonella for Powders, the method for Primary Production samples

SF: Salmonella for Feces, the method for fecal samples

Catalogue numbers

Salmonella Velox	SV96, SV960
Salmonella Velox 4 Powder	SP384
Salmonella Velox 4 Feces	SF384
Salmonella Velox SE+ST	SV96SEST SV960SEST
Salmonella Velox 4 Powder SE+ST	SP384SEST
Salmonella Velox 4 Feces SE+ST	SF384SEST

Other kits in this field

• Listeria Velox

Kit for 384 lysis and qPCR reactions for multiplex detection of Listeria spp. and/or Listeria Monocytogenes

• Salmonella Velox Cocoa

Kit for 96/960 lysis and qPCR 384 reactions for multiplex detection of Salmonella spp. and/or S. Enteritidis +S. Typhimurium (SE+ST)

For more info visit www.dna-diagnostic.com or contact Tel +45 87 32 30 50 · info@dna-diagnostic.com

DNA DIAGNOSTIC

DNA DIAGNOSTIC A/S · VOLDBJERGVEJ 14 · 8240 RISSKOV · DENMARK

About us DNA Diagnostic A/S is a Danish biotech company established in 1992. DNA Diagnostic develops and manufactures qPCR test kits for rapid identification of pathogenic microorganisms. DNA Diagnostic also makes CE-IVD kits for detecting leukemia related translocations. DNA Diagnostic is ISO 13485 certified.